


# Teruflex<sup>®</sup>

Multi-Function Hand Stripper  
& Aluminum Sealing Clips



## Multi-function Hand Stripper & Aluminum Sealing Clips

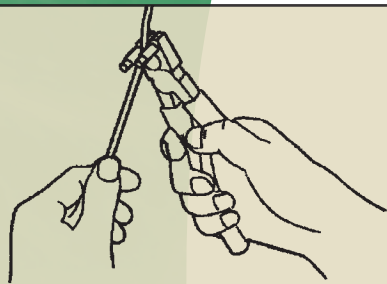


### Stripping

From needle end, strip toward bag.  
To strip, pull tube or push rollers.




### Sealing

Place sealing clip onto folded tubing.

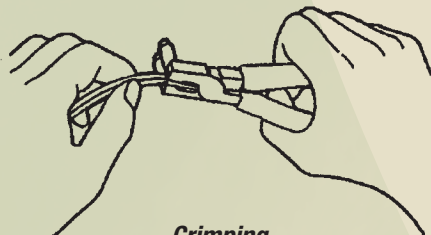

### Crimping

Crimp firmly until jaws close.

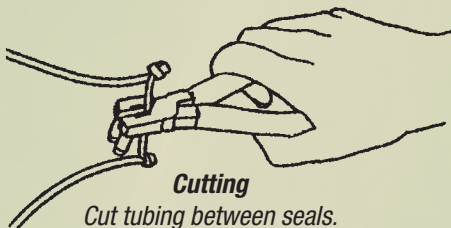

### Cutting

Cut tubing between seals.

### TERUFLEX® Multi-function Hand Stripper

- Does the work of three – strips, cuts and seals
- Spring-loaded design means no need to pull handles open after each strip, cut or seal
- Easy-to-grip, durable handle for secure control

### Aluminum Sealing Clip

- Ring design for easy placement onto folded tubing
- Made of aluminum for easy, secure sealing

Product Code	Description	Packaging
1ME-AC133	Hand Sealer	1 unit/case
1ME*ACS121	Aluminum Sealing Clip	1,000 units/case

### TERUMO MEDICAL CORPORATION

2101 Cottontail Lane  
Somerset, New Jersey 08873  
Tel: 732-302-4900  
Toll Free: 800-283-7866  
Fax: 732-302-4965  
Fax orders to: 800-411-5870  
[www.terumotransfusion.com](http://www.terumotransfusion.com)

TERUFLEX® and © TERUMO are registered trademarks of Terumo Corporation, Tokyo, Japan.  
All rights reserved. Copyright 2008 Terumo Medical Corporation.

1/08