

# SUCTION CANISTERS, LINERS & ACCESSORIES



# QUICK-FIT™ LINER SYSTEM



Now with Built-in Stopcocks

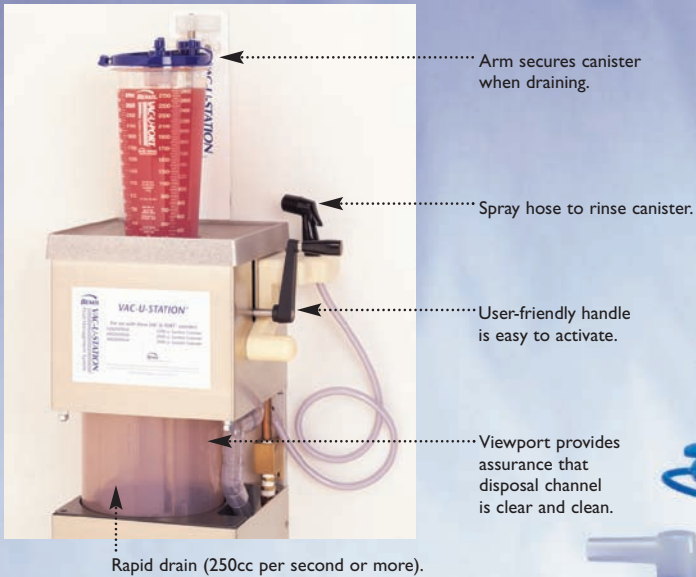
**QUICK-FIT LINER SYSTEM** provides easy insertion and safer removal of liners.

	DESCRIPTION	CATALOG NO.	PACKAGING	LINER SYSTEM
	<b>1500cc Quick-Fit Suction Liner</b> with attached lid	1504	50/case	 <p>Preassembled, sealed lid saves time and eliminates confusion.</p> <p>Handle offers easy removal and disposal.</p> <p>Innovative ribbed design and superior tear strength allow easy insertion and removal.</p> <p>Color coded lids and brackets allow easy visual selection.</p> <p>Liner unit comes compressed for packaging, for maximum use of storage space.</p>
	<b>3000cc Quick-Fit Suction Liner</b> with attached lid	3004	50/case	
	<b>1500cc Reusable Outer Canister</b> with attached bracket	1500 10 1500 01	10/case 1/case	
	<b>1500cc Reusable Outer Canister</b> with attached stopcock bracket	1500SC 10 1500SC 01	10/case 1/case	
	<b>3000cc Reusable Outer Canister</b> with attached bracket	3000 10 3000 01	10/case 1/case	
	<b>3000cc Reusable Outer Canister</b> with attached stopcock bracket	3000SC 01	1/case	



# VAC-U-PORT™ CANISTER SYSTEM

## Vac-U-Station™ Infectious Waste Management System


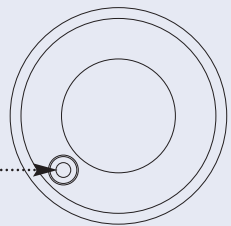




Secure disposal into a standard sanitary sewer system.

Stepped patient port on Vac-U-Port, Hi-Flow and Quick-Fit lids allows direct patient tubing connection.




**VAC-U-PORT CANISTER SYSTEM** directs suction canister contents directly into a sanitary sewer system.

DESCRIPTION	CATALOG NO.	PACKAGING	BOTTOM VIEW OF THE VAC-U-PORT CANISTER
 <b>1200cc Vac-U-Port Suction Canister</b> with Aerostat® filter, float valve shutoff	1202DSP 050	48/case	 <p>The canister has a solid base until it's ready to be drained. This patented drainage port is opened by the Vac-U-Station unit to form a pathway for the drainage process.</p>
 <b>2000cc Vac-U-Port Suction Canister</b> with Aerostat® filter, float valve shutoff	2002DSP 050	48/case	
 <b>3000cc Vac-U-Port Suction Canister</b> with Aerostat® filter, float valve shutoff	3002DSP 050	36/case	
<b>Vac-U-Station™ Unit:</b> (47.1" H x 18.7" W x 10.8" D) (See photo above)	VAC2000	1/case	

# HI-FLOW RIGID CANISTERS



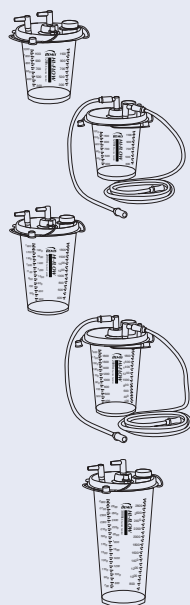**HI-FLOW RIGID CANISTERS** designed to meet high volume collection and disposal needs.



## DESCRIPTION

**800cc Hi-Flow Canister**  
with Aerostat® filter, float valve shutoff

**1200cc Hi-Flow Canister**  
with Aerostat® filter, float valve shutoff

**1200cc Hi-Flow Canister**  
with 18" (46cm) and 6' (2m)  
connecting tubes

**2000cc Hi-Flow Canister**  
with Aerostat® filter, float valve shutoff

**2000cc Hi-Flow Canister**  
with 18" (46cm) and 6' (2m)  
connecting tubes

**3000cc Hi-Flow Canister**  
with Aerostat® filter, float valve shutoff

**3000cc Hi-Flow Canister**  
with 36 count pkg. of 18" (46cm) connecting tubes

## CATALOG NO.

8002 055

484410

485410

494410

495410

3002 055

3025 055

## PACKAGING

70/case

48/case

48/case

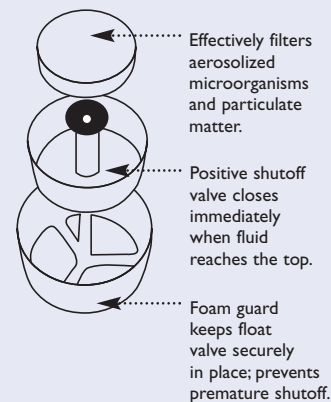
48/case

48/case

36/case

36/case

## AEROSTAT® HI-FLOW BACTERIAL FILTER





# HYDROPHOBIC RIGID CANISTERS

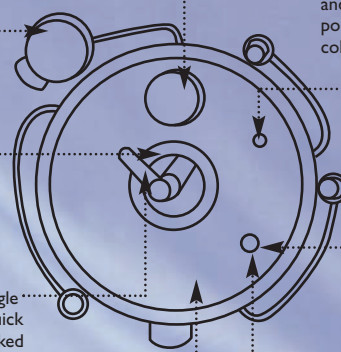
## Canister Lid Features



Attached port caps eliminate cap loss and allow quick, secure capping.

Vacuum port provides optimum flow rate and easy hookup.

Translucent right angle connectors allow quick hookup, prevent kinked tubing and allow for visual inspection.



Self-sealing lid provides a leak-free seal when vacuum is applied.

Large pour spout allows easy emptying and provides an adaptor port for the specimen collection accessory.

Tandem port for multiple canister connections.

Stepped patient port on Vac-U-Port, Hi-Flow and Quick-Fit lids allows direct patient tubing connection.

Large patient port ensures greater suction for maximum site clearance.

Hydrophobic 800cc canister lid features self-sealing lid, pour spout with attached cap, vacuum port and patient port. Right angle connectors are not required.

**HYDROPHOBIC RIGID CANISTERS** provide superior protection and maximum performance.

### DESCRIPTION

### CATALOG NO.

### PACKAGING

### HYDROPHOBIC SHUT-OFF FILTER



**800cc Hydrophobic Canister**  
with hydrophobic shutoff filter,  
built-in critical measure

424410

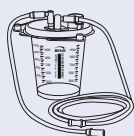
100/case



**1200cc Hydrophobic Canister**  
with hydrophobic shutoff filter

434410

48/case



**1200cc Hydrophobic Canister**  
with 18" (46cm) and 6' (2m)  
connecting tubes

438410

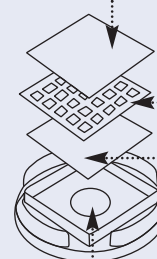
48/case



**2000cc Hydrophobic Canister**  
with hydrophobic shutoff filter

435410

48/case



Hydrophobic filter provides final filtration and overflow protection.

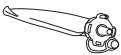

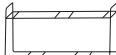
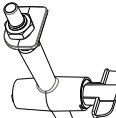
Breathing screen separates filter stages for improved suction.

Prefilter protects hydrophobic filter from erroneous shutoff.

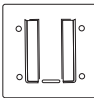
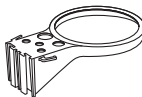
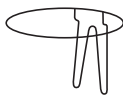








Tapered protective housing guards shutoff filter from splash or foam.

# ACCESSORIES









## ACCESSORIES

DESCRIPTION	CATALOG NO.	PACKAGING
 Specimen collection adaptor	533810	40/case
 Critical measurement insert for 1200cc, 2000cc and 3000cc canisters	534010	20/case
 Transport Tray	530000	10/case
 Stopcock for Rigid Canisters	532010	4/case






## BRACKETS

 Wall plate	530510	12/case
 Plastic ring for wall plate for 1200cc, 2000cc and 3000cc canisters	532710	12/case
 Wire ring for wall plate for 1200cc, 2000cc and 3000cc canisters	532810	12/case
 Wire ring for wall plate for 800cc canister	534510	12/case
 Low-boy metal ring for wall plate for 1200cc, 2000cc and 3000cc canisters	533010	12/case
 Wire ring for regulator mount for 1200cc, 2000cc and 3000cc canisters	532910	12/case
 Wire ring for regulator mount for 800cc canister	534610	12/case
 Metal regulator bracket/handle for roll dolly stands	530310	1/case
 Quick-connect regulator mount for 800cc hydrophobic canister	534410	1/case
 Regulator bracket for Quick-Fit liner	538510	12/case
 Wire hanger for Quick-Fit liner	539010	12/case

## TANDEM HOOKUPS

DESCRIPTION	CATALOG NO.	PACKAGING
 Single liner stopcock manifold	538110	25/case
 Two-liner stopcock manifold (not shown)	538210	5/case
 Two-canister stopcock manifold	535510	5/case
 Four-liner stopcock manifold	538410	5/case
 Four-canister stopcock manifold	535010	5/case
 Four-liner manifold - no stopcocks (not shown)	539410	5/case
 18" nonsterile connecting tube	536510	100/case
 6' nonsterile connecting tube	536010	25/case

## ADAPTING HARDWARE

 Male 1/8" N.P.T. x Tubing Barb	530910	1/case
 Male 1/8" N.P.T. x DISS	531910	1/case
 Male 1/8" N.P.T. x DISS x Universal Tubing Barb	531210	1/case
 Female DISS x Tubing Barb	531110	1/case
 Female DISS x Universal Tubing Barb	531410	1/case

## STANDS

A	Four-canister Hi-low stand	533910	1/case
	Four-canister stand for liners	534910	1/case
	Eight-canister adaptor (not shown)	538010	1/case
B	Single roll dolly	532410	1/case
C	Double roll dolly	532310	1/case

